

# Denyer XRT/1500



## Operating Table

DENYERS

# The Denyer *XRT1500*



Prone / Supine



Beach chair



Trendelenburg



Reverse Trendelenburg



Lithotomy



Low Lithotomy



Lateral, Kidney etc



Lateral Tilt

# Surgical Operating table

## Introduction:

The new XRT1500 series continues the Denyer tradition of high quality, reliable and well priced surgical tables which incorporate the experience and feedback gained in over 120 years of working directly and closely with health professionals.

## Surgery types:

The XRT1500 is designed for basic surgical procedures such as Gynaecological, Urological, Orthopaedic, Plastic, Ophthalmic and General, and is particularly well suited to day procedure hospitals.

## Controls, ease of use and cleaning:

The table is primarily controlled by a lightweight, easily read hand piece which controls height, Trendelenburg, lateral tilt, back section up/down and foot brakes. Return-to-level is accomplished with a single button which brings all powered sections back to level with a single press. The head and foot sections are interchangeable and manually actuated and the table can be ordered with either a single or split leg plates.

Ease of use and practicality were driving issues in design and to this end the hand control has simple, very legible buttons and all controls have large articulation ranges. All main surfaces are high quality stainless steel with a high gloss phenolic resin top surface which makes for easy and efficient cleaning.

## Power supply:

The powered movements are actuated by an hydraulic pump driven by rechargeable batteries. The batteries have enough capacity for several days of operations, but should the batteries ever go flat the table will operate directly from mains power.

## Radiolucency:

The table has been designed to accommodate the latest surgical techniques, allowing exceptional Surgeon access and high levels of radiolucency to all parts of the body with image intensifiers.

## Weight limits:


With obesity becoming such a major issue, the XRT1500 is designed to accommodate patients of up to 220kg when centrally positioned, or up to 180kg in lithotomy or any other positions where the patient is significantly offset from the centre of the table.

## Accessories:

A large range of accessories, including lightweight extensions, positioning accessories, stirrups, clamps and many others is available, and Denyers also have a full range of specialist operating tables – please contact your local Denyer Representative or visit our web site at [www.denyers.com.au](http://www.denyers.com.au) for more information

<b>Specifications:</b>		
Height	660mm (26") min	1020mm (40") max
Width of top	523mm (21")	575mm (23") incl. rails
Length of top	2055mm (81")	
Trendelenburg	30°	30° reverse Trend
Flex / reflex	90°	220°
Lateral tilt	20°	
Back up / down	90°	40°
Leg up / down	30°	95°
Head up / down	45°	90°
Patient weight	220kg (485lbs) max	180 kg (400lbs) lithotomy
Because of our policy of continual product development, all specifications are subject to change without notice.		

## History:

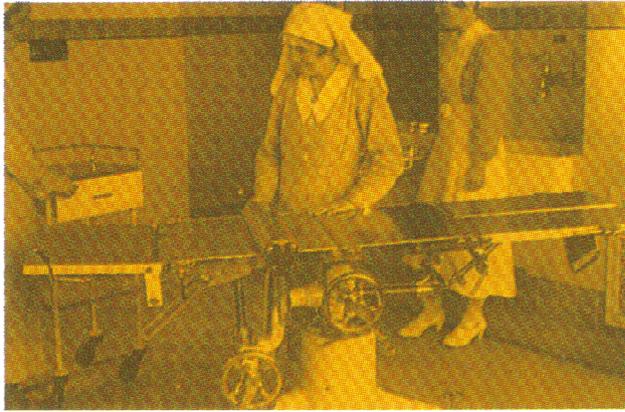

Since its inception in 1889, Denyers has always followed the same philosophy—to produce the highest quality product range without compromise and to be leaders in their design and manufacture. The company strives to improve its product range to the benefit of end users, adjusting to the requirements of the international market where

differences in surgical techniques place other demands on equipment design.

Denyers have an historical commitment to excellence and product life, which in turn ensures lifetime protection for the customer's investment.

With 120 years of service to the medical profession, Denyers have always listened to their customers' wants and needs and this is why Denyers have become a premier manufacturer of surgical operating tables.



Today Denyers is small enough to be flexible yet large enough to implement change quickly with a constant eye on the user's benefits and an unwavering policy of engineering excellence.

**DENYERS**  
INTERNATIONAL

Email [sales@denyers.com.au](mailto:sales@denyers.com.au)  
Web [www.denyers.com.au](http://www.denyers.com.au)

Melbourne	03 9555 5722
Sydney	02 9627 3299
Brisbane	07 3390 7011
Adelaide	08 8354 0933
Perth	08 9284 1033
International	+61 3 9555 5722

42—46 Nelson Street  
Moorabbin Vic  
Australia 3189



Your local Representative:

