

# XRT4200 Specialist Table



EXPERIENCE THE SUPERIOR QUALITY, FLEXIBILITY & EXCEPTIONAL IMAGING

[WWW.DENYERS.COM.AU](http://WWW.DENYERS.COM.AU)

## Gynaecology, Urology



PROUDLY MADE IN AUSTRALIA

# **DENYERS**

INTERNATIONAL

# our solution XRT4200

The new XRT4200 table continues the Denyer tradition of high quality, reliable and well priced surgical tables which incorporate the experience and feedback gained over 130 years of working with Surgeons, Anaesthetists, Nurses and other medical specialists. The XRT4200 is Australian designed and assembled.

## XRT4200 KEY FEATURES

- Specially designed for gynae work
- Specifically designed to allow transfer to table directly into lithotomy position - no repositioning of patient after transfer
- Anaesthetist retains articulation control of head throughout entire time patient is on the operating table
- Equally capable as a general operating table, with excellent surgeon and C-arm access in prone, supine, lateral



\* The lightweight leg section & urology section are standard with XRT4200 Gynae/Uro table configuration

\* If the table is required for General use a standard leg section & seat section will be required.

## SURGERY TYPES

The XRT4200 is suitable for all disciplines of surgery, particularly Gynae/Urology with interchangeable lower trunk, leg and head sections allowing customisation for the surgeon yet still providing full access for anaesthetists.

## RADIOLUCENCY


The table has been designed to accommodate the latest surgical techniques, with a back-lit hand piece to assist in dimly lit endoscopic procedures and with exceptional radiolucency for ease of access to all parts of the body with image intensifiers.

## CONTROLS, EASE OF USE AND CLEANING

The table is primarily controlled by a lightweight, easily read hand piece which controls height, Trendelenburg, lateral tilt, back section up/down and foot brakes. An optional foot control is available for control of height, Trendelenburg and lateral tilt.

All main surfaces are high quality stainless steel with a high gloss phenolic resin top surface which makes for easy and efficient cleaning.

## XRT4200 | SPECIFICATIONS

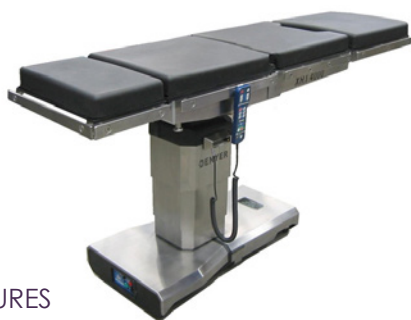


### KEY FEATURES

- Substantially greater patient weight capacity
- Great C-arm access for Gynae, Urology
- Infra red remote option, Optional lower height\*
- Australian designed and manufactured

Height	660mm (min)	1120mm (max)
*Low Height Version	*620mm (min)	*1040mm (max)
Width of top	523mm	575mm (incl. rails)
Length of top	2055mm	
Trendelenburg	30°	30° reverse
Flex / reflex	220°	90°
Lateral tilt	20° Left / Right	
Back up / down	90° up	40° down
Leg up / down	30° up	95° down
Head up / down	45° up	90° down
Patient weight	500kg max	280kg lithotomy



## GYNAE/URO TABLE CONFIGURATION

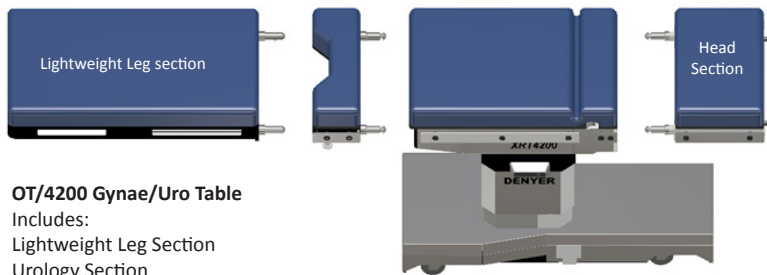

### OT/4200 Gynae/Uro Table

Includes:  
Lightweight Leg Section  
Urology Section  
Head Section

## GENERAL USE TABLE CONFIGURATION

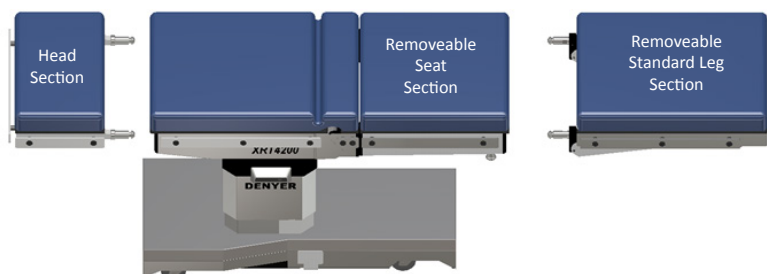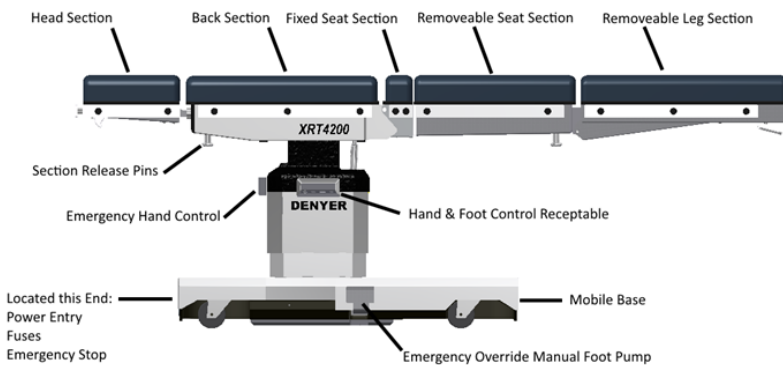

## OPTIONAL TABLE EXTENSIONS



## TABLE ACCESSORIES



**AM/O-UC4 Uro Catcher**  
AM/O-UC31 Sterile Bags  
AM/O-UC32 Non-Sterile Bags



Trend Lateral Height  
**OTA/4001 Surgeon's Foot Control**



**AM/A-60001 Hug-U-Vac**  
Available as Standard or Steep Trend in Standard and large sizes.



**AM/A-63500 Shoulder Braces**  
Support without pressure on Brachial Plexus. Sold as a pair.  
\* Requires 2 x rail clamps

## STIRRUPS

**AM/O-YFASI Yellofin Stirrups**  
Gas assisted lithotomy Stirrups  
Patient weight limit 159kg

**Requires:**  
AM/O-AMA68S0.6S  
Reuseable boot pads

2 x AM/O-RC2 rail clamps



**AM/O-YFESP Pinkfin Stirrups**  
Gas assisted lithotomy Stirrups  
Patient weight limit 227kg

**Requires:**  
AM/O-AMA68S0.6S  
Reuseable boot pads

2 x A M/O-RC2 rail clamps



**AM/O-UFAS Ultrafin Stirrups**  
Gas assisted lithotomy Stirrups  
Patient weight limit 363kg

Reusable pads included

2 x A M/O-RC2 rail clamps





# DENYERS

INTERNATIONAL



## EXCELLENCE SINCE 1889

Since its inception in 1889, Denyers has always followed the same philosophy—to produce the highest quality product range without compromise and to be leaders in their design and manufacture.

Denyers have a historical commitment to excellence and product life, which in turn ensures lifetime protection for the customer's investment. With 130 years of service to the medical profession, Denyers have always listened to their customers' wants and needs and this is why Denyers have become the leader in our chosen field.

Today Denyers is small enough to be flexible yet large enough to implement change quickly with a constant eye on the user's benefits and an unwavering policy of engineering excellence.

Denyers are certified to ISO9001 & ISO13485 by SGS

DEN\_002

DENYERS INTERNATIONAL LTD

Freecall: 1800 DENYER (1800 336 937)

Head Office

42 - 46 Nelson Street, Moorabbin, Victoria 3189, Australia

Phone: +61 3 9555 5722

Fax: +61 3 9555 4722

Email: [sales@denyers.com.au](mailto:sales@denyers.com.au)

[www.denyers.com.au](http://www.denyers.com.au)

DENYERS NZ LTD

Tollfree: 0800 DENYERS (0800 336 9377)

Phone: +64 9 215 6671

Email: [sales@denyers.co.nz](mailto:sales@denyers.co.nz)

[www.denyers.co.nz](http://www.denyers.co.nz)